

veloce

Intense Pulsed dual Light System for Multiple Indications

Skin Tightening with Advanced Pulsed Light - N-Tight Mode

Consistent High Peak Power, Powered by Dual Lamps

Improved Cooling System for Safety and Patient Comfort

Broad Spectrum Light with 10 Interchangeable Filters Covering Multiple Treatment Indications



ParagonCare
Aesthetics

One Platform. Multiple Indications. Smarter ROI.

IPDL(Intense Pulsed Dual Light) System with Broad Band Light

VELOCE is an advanced IPDL system designed for Skin tightening, pigmentation, hair removal, acne and skin rejuvenation. It delivers stable and accurate energy output throughout treatment, ensuring consistent clinical outcomes. With broad treatment versatility, it effectively supports a wide range of indications across different skin concerns helping clinics' ROI improvement.


Veloce's Core Technology

Skin Tightening with Advanced Pulsed Light - N-Tight Mode

Broadband near-infrared light delivers controlled thermal energy to the dermis, supporting collagen contraction and remodeling. Combined IPDL and N-Tight technology enhances skin firmness while addressing tightening.

Consistent High Peak Power, Powered by Dual Lamps

While some IPL systems may show energy variability during continuous operation, the dual-lamp architecture supports stable, consistent delivery at 20 sub-pulses per second, generating 320 sub-pulses within 16 seconds (1,080 J/sec) and covering up to a 45 × 15 mm spot size for efficient, high-throughput treatments.

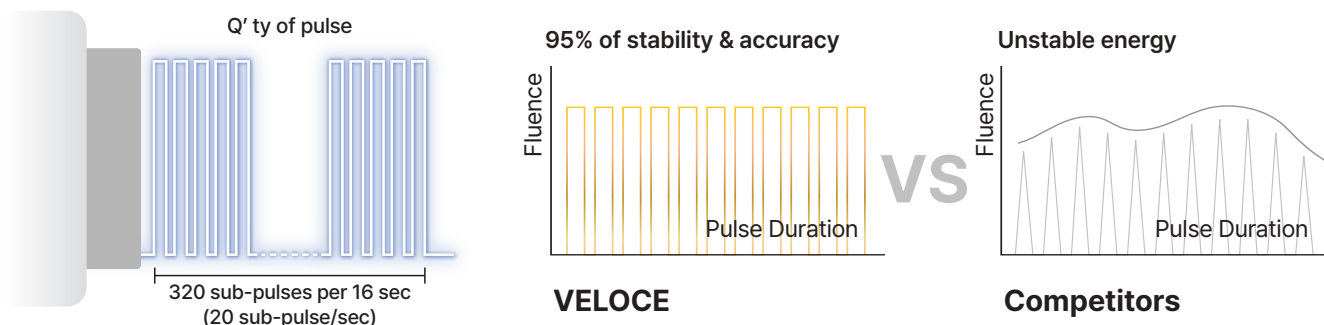


Improved Cooling System for Safety and Patient Comfort

The sapphire tip's powerful cooling ability reduces the skin temperature to almost 0°C within 3 minutes. The rapid and impressive cooling control protects the area to be treated from the damage and keep patients safe during the entire treatment.



Broad Spectrum Light (420 -1200 nm) for Multiple Indications

With 10 interchangeable filters, VELOCE enables customizable treatments for diverse skin concerns with greater efficiency.

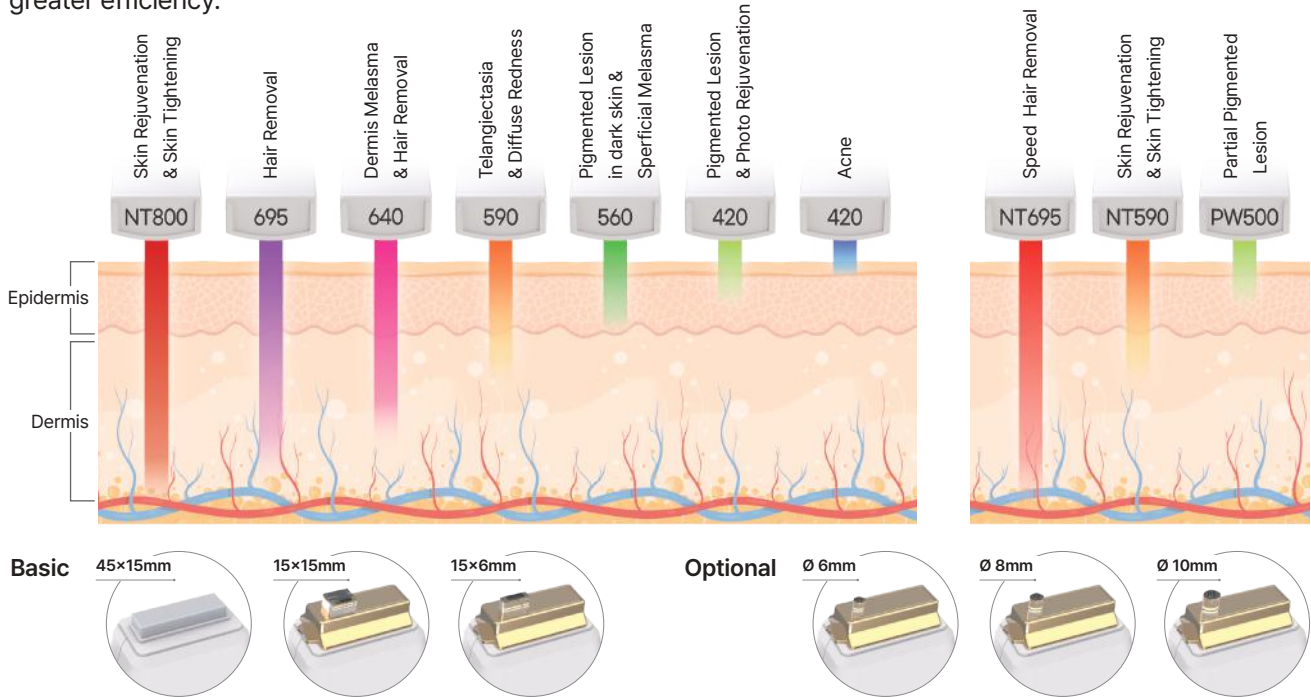

Clinical Results

Neck Tightening
(VELOCE + Mesotherapy) (Courtesy of Dr. Yelda Kapicioğlu)



Baseline



After 1 session

Sun-Damaged Scars & Freckles
(Courtesy of Dr. Erol Koç)



Baseline



After 1 session

Nevus Flammeus
(Courtesy of Dr. Yelda Kapicioğlu)



Baseline



After 1 session

Rosacea & Telangiectasias
(Courtesy of Dr. Erol Koç)



Baseline



After 1 session

veloce

Specification

Wavelength	420 ~ 1200nm
Filters	420nm, 515nm, 560nm, 590nm, 640nm, 695nm, NT800nm, PW500nm (Option) NT695nm, NT590nm
Fluence	30J/cm ² Max.
Max Energy	202J
Pulse Width	5 ~ 200ms
Repetition Rate	Max. 1Hz
Cooling Temperature	0 ~ 25°C (Room Temp 25°C)
Spot Size	45 × 15mm, 15 × 15mm, 15 × 6mm (Option) Ø6mm, Ø8mm, Ø10mm
Weight	85kg
Dimensions	480(W) x 826(D) x 1094(H)mm
Power Consumption	220VAC, 50/60Hz, 4.3KVA



ParagonCare
Aesthetics

Classys Distribution Partner of Australia & New Zealand



Manufacturer Classys Inc. CLASSYS 208 Teheran-ro
Gangnam-gu, Seoul, Republic of Korea, Seoul KOR.
Australian Sponsor Classys Australia Pty Ltd

Book a demo now

 **1300 633 603**

 **aesthetics@paragoncare.com.au**

   **@paragoncare.aesthetics**